

Reflex XL & Reflex XXL shower toilet chair


Instruction card

01. Safety recommendations

- Read instruction for use before use (www.lopital.com).
- Maximum load 227 kg / 500.45 lbs (Reflex XL); 275 kg / 606.3 lbs (Reflex XXL).
- Make sure that the chair is on the brake during client transfers from and to the device.
- Place both feet on the foot supports.
- Place the arms on the armrests.
- Only relocate and move the chair in the lower position.




02. Control panel

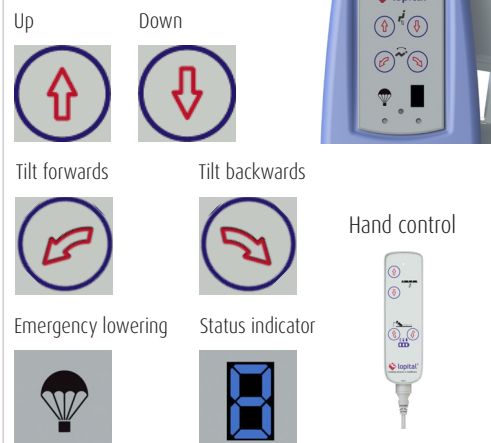


03. Armrest

1. Release
2. Press the blue button
3. Swing upwards

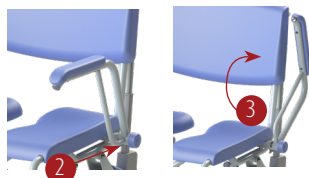

5100 5900
Reflex XXL

5100 5950
Reflex XL

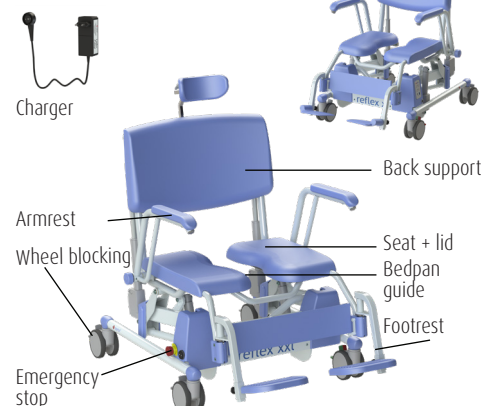

04. Use

1. Remove the Reflex XL/Reflex XXL from the charger, if necessary.
2. Position the chair correctly and, if necessary, fold the footrests up or remove the footrests from the Reflex XL/Reflex XXL.
3. Apply the brakes to the chair!
4. The client can now take a seat in the chair. Assist the client as needed. If desired, the chair can be tilted forward by 10° to facilitate seating.
5. Place both feet on the footrests with the calves against the calf strap.


05. Footrest

Both footrests must be swung away before the client sits down in the chair. Ensure that the client's feet are resting on both footplates while the chair is being moved.

The footrests of the Reflex XL/XXL are height-adjustable, and the footrests of the Reflex XXL are also width-adjustable.

06. Charging batteries

After using the Reflex XL/XXL, connect it to the charger to maintain the batteries.

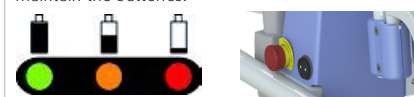

Off
Green
Orange
Red

Sleep setting
Chair ready to use
Charging
Battery is empty


The Reflex XL/XXL will not work during charging.

07. Toilet function

- Remove the lid from the seat before the client with a disability takes a seat in the chair.
- *In combination with toilet:* the Reflex XL/XXL can be moved over almost every toilet. Move the chair upwards or downwards to achieve the right position.
- *With bedpan:* the Reflex XL/XXL can also be used as a toilet chair. Simply slide the stainless steel / plastic bedpan under the chair using the guides.


08. Emergency stop & lowering button

1. Activate the emergency stop in the event of an emergency where the upward or downward movement must be stopped immediately.
2. When the lowering button is activated in an emergency, the chair lowers to its lowest position.

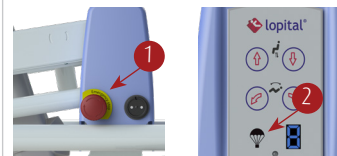

09. Wheel blocking

1. Wheel brake.
2. Unlocking.
3. Directional lock.

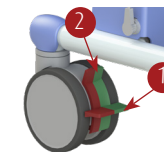

10. Cleaning and maintenance

The Reflex XL/XXL must be cleaned after every use, with water and regular household cleaning agents. Do not use any abrasive or other aggressive cleaning products. If necessary, the PUR components can be disinfected with a solution of 70% ethanol or propanol. The complete maintenance schedule and cleaning protocol can be found in the instructions for use (www.lopital.com).

According to the regulation MDR (EU) 2017/745 for medical devices, Reflex XL/XXL must be maintained in line with the manufacturer's instructions. Those instructions can be found in the maintenance protocol, see website Lopital (www.lopital.com).

Important!

This quick instruction for use may not be considered a replacement for the instruction for use!

In the event of a breakdown, contact Lopital.
Tel: +31 (0) 13-5239300
E-mail: info@lopital.nl

Instruction for use
For the full instruction for use go to: www.lopital.com

