



Maintenance protocol

**Tango XL & Tango XXL
shower toilet chair**



Tango XL shower toilet chair: 5100 5750



Tango XXL shower toilet chair: 5100 5850

Important:

The Tango XL / Tango XXL must be serviced and inspected at least once a year.

Service and maintenance of the Tango XL / Tango XXL may only be carried out by persons qualified by Lopital and using original spare parts.

Modifications to the construction may compromise the safety of the Tango XL / Tango XXL. Furthermore, this will void Lopital's liability and warranty conditions. As a result, the Tango XL / Tango XXL will no longer comply with the MDR (EU) 2017/745 regulation for medical devices.

Any inspection or repair must be recorded in the logbook.

The average technical service life of the Tango XL / Tango XXL is 10 years, provided that daily and periodic maintenance is carried out.

Make sure that all screws, bolts, and nuts are checked and tightened if necessary.

Manufacturer

Lopital Nederland B.V.

Laarakkerweg 9, 5061 JR Oisterwijk

T +31 (0)13 5239300

www.lopital.com, info@lopital.nl

MAIN5850-2026-01_EN

Table of contents

1.	Overall condition of the chair.....	4
2.	Wheels	5
3.	Armrests.....	6
4.	Footrest.....	7

1.	Overall condition of the chair
2.	Wheels
3.	Armrests
4.	Footrest

C = Check
 S = Lubricate
 A = Recommended replacement (or replace parts)
 V = Replace (or replace parts)
 R = Clean



Month	12	24	36	48	60	72	84	96	108	120
Overall condition	C	C	C	C	C	C	C	C	C	C
PUR	C	C	C	C	A	C	C	C	C	A
Wheels	C	C	C	C	C	C	C	C	C	C
Armrests	C	C	C	C	C	C	C	C	C	C
Footrests	C	C	C	C	C	C	C	C	C	C

1. Overall condition of the chair

Month	12	24	36	48	60	72	84	96	108	120
	C	C	C	C	A	C	C	C	C	A

Required tools

- Paint pen A77A2

- Check that all stickers are still present and legible.
 - Serial number sticker:
 - Is the serial number still legible?
 - Is the data matrix still legible for a scanner?
 - Is the CE logo still visible?
 - Inspection sticker:
 - Is the inspection date still legible?
 - QR sticker:
 - Is the QR code still legible for a QR scanner?
- Visually inspect all painted parts for damage and rust. Touch up with paint if necessary.
- Check the PUR parts for severe discolouration (especially brown) and damage.
- Check that all screws are present and tight. This includes the guide!



Figure 1 Stickers



For hygiene reasons, it is strongly recommended that the seat and backrest be replaced every 5 years.

Included (can be removed without tools):

- 2 seat halves
- Footrest (not a rejection point)
- Calf strap
- Tango XXL backrest

2. Wheels

Month	12	24	36	48	60	72	84	96	108	120
	C	C	C	C	C	C	C	C	C	C

Required tools
- Allen key 8 mm

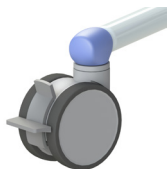


Figure 2 Tango XXL wheel



Figure 3 Tango XXL wheel exploded view



Figure 4 Tango XL wheel

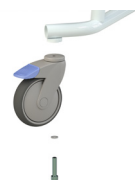


Figure 5 Tango XL wheel exploded view

- Check the wheel attachment (play).
- Check the free wheel + brakes.
- Check that all 4 wheels are touching the ground. If not, check that the chair tracks straight.
- Check that the brake tab is intact.
- Check the tread.



If any of these conditions are not met, replace the wheel.

3. Armrests

Month	12	24	36	48	60	72	84	96	108	120
	C	C	C	C	C	C	C	C	C	C

Required tools
– Allen key 4 mm

- Check the operation of the armrests.
- Check the armrest for damage, rust and clearance.
- Check that all screws are present and securely tightened.
- Check the width adjustment of the armrests on the Tango XXL.

For an exploded view, visit www.lopital.nl



Figure 6 Tango XL armrest



Figure 7 Tango XXL armrest

4. Footrest

Month	12	24	36	48	60	72	84	96	108	120
	C	C	C	C	C	C	C	C	C	C

Required tools

- Allen key 4 mm

- Check the footrest for damage.
- Only Tango XL: fold up the footplates (they should remain in the folded position without falling back automatically).
- Only Tango XXL: check the width adjustment (swing-away function) of the footrest.
- Check that the footplates are level and slightly angled upwards towards the centre. If not, also check the frame.
- Check that the footrest is securely fastened at different heights. Try various height settings.
- Is the pull knob present?
- Only Tango XL: is the footplate cap present?



Footplate up

Figure 8 Tango XL footrest



Figure 9 Tango XXL footrest

For questions about the maintenance protocole, please contact Lopital (telefoon +31 (0)13 5239300 or email info@lopital.nl). You must provide Lopital with the serial number (red border). See the serial number sticker on the product for this.



lopital®
 Laarakkerweg 9
 5061 JR Oisterwijk
 The Netherlands
 T: (031)13-5239300

IPX-
 Weight: -- kg
 Type battery: ---
 Drives operates: -- V
 Charge: ---V--- mA
 Non continuous: -----

MD -----

REF -----

SN -----

MAX
LOAD ----- KG/---LBS

 -----



(01)-----
 (11)-----
 (21)-----

CE

Revision -



If caregivers have a choice



For more information visit our website: www.lopital.nl

Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk, Postbus 56,
5060 AB Oisterwijk, Tel +31 (0)13.52.39.300, E-mail info@lopital.nl,
www.lopital.nl

For more information visit our website: www.lopital.be

Lopital België B.V., Antwerpsesteenweg 124, 2630 Aartselaar
Tel +32 (0)3/870.51.60, E-mail info@lopital.be, www.lopital.be

Copyright ©2026 Lopital. All rights reserved.



Green Deal



**MADE IN
THE NETHERLANDS**



ISO CERTIFIED
9001 - 13485 - 14001



**Scan for
information**