

Reflex

Tilted for comfort, built for reliability

5100 5600

Reflex 5600 shower toilet chair with height-adjustable footrest

This electric, height-adjustable shower toilet chair is designed to provide safe, dignified care with maximum comfort for both the user and caregiver. The integrated seat mimics the support and feel of a regular chair while offering essential functionality such as electric height control and up to 30° tilt. This makes it ideal for both washing and caring activities in a range of care environments. The multifunctionality of the design ensures that personal hygiene care can be delivered ergonomically and efficiently.

For additional versatility, the chair includes detachable and foldable footrests with calf protection, allowing for easier transfers and improved leg support. Its foldable armrests facilitate lateral movement, while the robust base frame is built for safe and stable positioning over a toilet or when using a bedpan or bucket. Smooth electric operation and intuitive controls make it a practical solution for caregivers and a comfortable, secure seat for users with mobility limitations.

Ergonomic Use. Equipped with smart handles and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.

Reliable & Long-Lasting. Built for daily care and backed by a 60-month warranty, offering a dependable performance.

Antibacterial & Corrosion-Resistant. The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

 Discover our tools in action through the Interactive Product videos on our website: www.lopital.com



Technical specifications

Max. load	135 kg / 300 lbs
Product weight	53 kg / 116.8 lbs
Seat width	47 cm / 18.5 inch
Seat depth	44 cm / 17.3 inch
Seat height (adjustable)	50 to 100 cm / 19.7 to 39.4 inch



Suitable for people with very slightly to very severely reduced mobility

5100 5600

Reflex 5600 shower toilet chair

Height and depth-adjustable head / neck support

Ergonomic handle for smooth and light operation

Ergonomic seating comfort

Guide rails for bedpan or toilet bucket

Strong and stable frame for maximum safety and support

Armrests with front closure for safety and security

Detachable and height-adjustable foot-rests and calf protection

Bumper to prevent damage to the frame

Charging point with magnetic coupling, no separate battery required.


Control panel.

Optional

7057 0508
Castors 125 mm
Directional.



6840 5661
Dual locking
armrests.

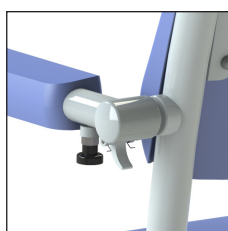

6950 5624
Platilon backrest.



Options

- Product no. 6830 5650 Straight armrests
- Product no. 6950 1104 Bedpan with lid plastic
- Product no. 6975 8044 Seatbelt 2-point
- Product no. 6975 8045 Seatbelt 4-point
- Product no. 6950 5624 Platilon backrest
- Product no. 6950 5625 Platilon seat U-model
- Product no. 7057 0508 Castor 125 mm Directional

Technical specifications

Total length base frame	106 cm / 41.7 inch
Total width base frame	70 cm / 27.6 inch
Turning radius	120 cm / 47.2 inch
Tilt setting	0 to 30 degrees
Electrical power	24V
Battery capacity	7Ah
Waterproofness	IPX4
Number of lifting actions	40 up-down movements with fully charged batteries
Head/neck support depth/height	7 - 12 cm / 2.8 - 4.7 inch

5100 5605

Reflex 5605 shower toilet chair

This height-adjustable shower toilet chair serves as a multifunctional solution for hygienic care, supporting clients with reduced mobility. It consists of a base frame and an ergonomically shaped seat that provides both comfort and functional stability. Gradual electric height and tilt adjustments allow caregivers to work in safe, ergonomic postures while enabling clients to remain supported and comfortable throughout care routines. The chair can also tilt backwards up to 30°, offering enhanced safety and relaxation.

Its design facilitates toileting and washing without unnecessary transfers, and a rail system that accommodates a bedpan or collection bucket. The sturdy frame is highly manoeuvrable and equipped with four braked castors for safety and stability during use. Whether used bedside, in a bathroom, or over a toilet, this device provides a seamless blend of convenience, hygiene, and well-being for both caregiver and client.



Head/neck support



Swivel castor with individual brake



Control panel.



Technical specifications

Max. load	135 kg / 300 lbs
Product weight	53 kg / 116.8 lbs
Seat width	47 cm / 18.5 inch
Seat depth	44 cm / 17.3 inch
Seat height (adjustable)	50 to 100 cm / 19.7 to 39.4 inch

Discover our tools in action through the Interactive Product videos on our website: www.lopital.com



Suitable for people with very slightly to very severely reduced mobility

5100 5605

Reflex 5605 shower toilet chair

Optional

6960 8026
Adjustable comfort leg support.

Ergonomic handles for smooth and light operation


Armrests with front closure for safety and security

Guide rails for bedpan or toilet bucket

Ergonomic seating comfort

Strong and stable frame for maximum safety and support

Bumper to prevent damage to the frame

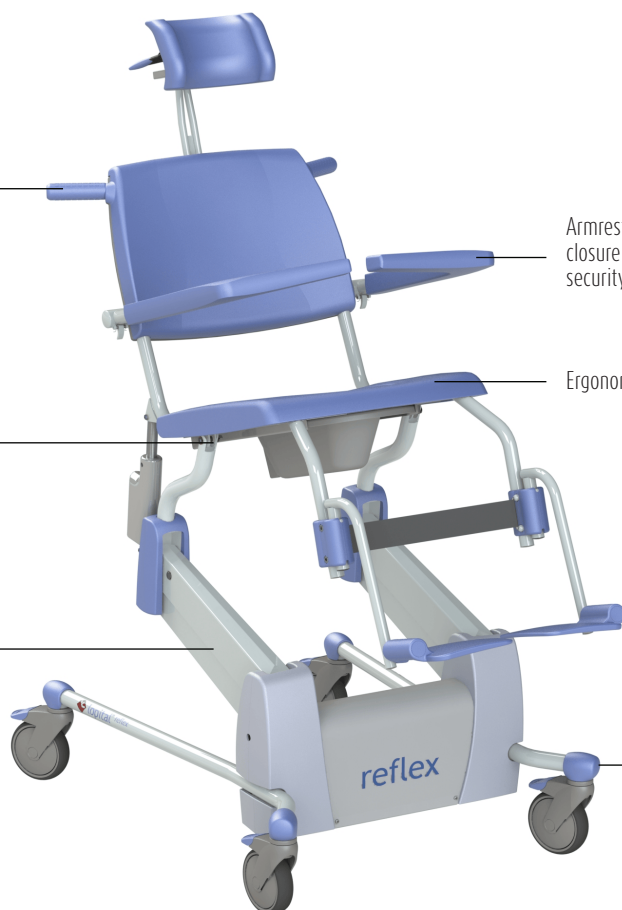

Technical specifications

Total length base frame	106 cm / 41.7 inch
Total width base frame	70 cm / 27.6 inch
Turning radius	120 cm / 47.2 inch
Tilt setting	0 to 30 degrees
Electrical power	24V
Battery capacity	7Ah
Waterproofness	IPX4
Number of lifting actions	40 up-down movements with fully charged batteries
Head/neck support depth/height	7 - 12 cm / 2.8 - 4.7 inch

Options

- Product no. 6830 5650 Straight armrests
- Product no. 6950 1104 Bedpan with lid plastic
- Product no. 6975 8044 Seatbelt 2-point
- Product no. 6975 8045 Seatbelt 4-point
- Product no. 6950 5624 Platilon backrest
- Product no. 6950 5625 Platilon seat U-model
- Product no. 6960 8026 Comfort leg support set
- Product no. 7057 0508 Castor 125 mm Directional

LFT5600-2025-01 EN



For more information also visit our website: www.lopital.com
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk
Tel +31 (0)13 52.39.300, E-mail info@lopital.nl, www.lopital.nl



ISO CERTIFIED
9001 - 13485 - 14001

MADE IN
THE NETHERLANDS



Scan for
information

All of the aids produced by Lopital comply with the applicable standards for Class 1 medical products (MDR (EU) 2017/745).

Copyright © 2025 Lopital. All rights reserved.