

# Marina

Ergonomic, durable and stylish in any healthcare environment

**Hydraulic: 6100 2010 / Electric: 6100 2060**  
**Marina Deluxe shower trolley**

The Marina shower trolley is designed for high-quality care in the lying position. This trolley supports caregivers in providing efficient, safe and dignified care, without physical strain. Thanks to the infinitely adjustable backrest, Anti-Trendelenburg position and height adjustment (electric or gas spring), the working posture is optimized – resulting in fewer physical complaints and more restful handling.

For clients, this means a stable, comfortable experience where safety and respect are paramount. The double-sided access makes care operations more flexible and reduces time per client, contributing to an efficient workflow.

The durable finish with antibacterial silver ion coating, corrosion-resistant frame and easy-to-clean surfaces provide a safe, hygienic solution that lasts. In addition, the sturdy castors with central brake and compatibility with patient elevators and glide rails ensure safe and smooth integration into any care environment.

Lopital Marina – comfort for the client, ergonomic ease for the caregiver.

**Standard built-in length Adjustment & Dual-sided Care Access.** Designed for 360° access with side rail controls on both sides and adaptability due to built-in length adjustment.

**Plug & Play Design.** Ready to use out of the box – intuitive, low-maintenance, and fast to implement in any care setting.

**Ergonomic Use.** Equipped with smart handles and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.

**Reliable & Long-Lasting.** Built to last with a 60-month warranty for long-term reliability and peace of mind.

**Antibacterial Coating & Corrosion-Resistant.** The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.



Made of 75% recycled aluminium



Lowest CO<sub>2</sub> footprint



Discover our tools in action through the Interactive Product videos on our website: [www.lopital.com](http://www.lopital.com)



Suitable for people with severely reduced to very severely reduced mobility

## Technical Specifications

High / low drive	Hydraulic or electric
Back adjustment	Gas spring or electric
Manual	Four-button hand control with hanging hook, IPX6
Control box	24VDC, max. 250 W, IPX6
Battery box	24V, 2.9 Ah, IPX6
Charger	100-240V AC, 50-60 Hz, max 1A, Wall charger

## Marina Deluxe




70205004  
Hand control  
with battery  
indication



70205003  
Battery box

### Optional

70200034 Pillow  
(a pillow is included with  
the basic models).

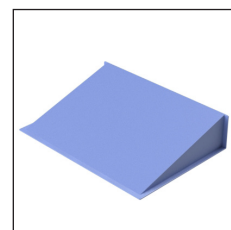

Art.no.	Art. description	Outer length	Outside width	Inner length	Inner width	High / low	Max. load	Product weight in kg
6100 2000	Marina Basic hydraulic shower trolley	186 206 226	81	170 190 210	70-80	45-95	227	76
6100 2010	Marina Deluxe hydraulic shower trolley	186 206 226	81	170 190 210	70-80	45-95	227	80
6100 2050	Marina Basic electric shower trolley	186 206 226	81	170 190 210	70-80	45-95	227	76
6100 2060	Marina Deluxe electric shower trolley	186 206 226	81	170 190 210	70-80	45-95	227	80

\*All sizes in this table are in centimetres.

Hydraulic: 6100 2000 / Electric: 6100 2050

## Marina Basic shower trolley

The Marina Basic is a functional and user-friendly shower trolley for recumbent care.

Unlike the Marina Deluxe, this basic model does not have a backrest. This model has four swivel castors with individual brakes, making it easier to use and maintain. The Marina Basic is available with hydraulic or electric height adjustment and is equipped with folding side rails for safe and comfortable care. An ideal solution for healthcare institutions looking for functionality, reliability and ease of use.



Hydraulic control pedal.



Swivel castor with individual brake.



Lowest CO<sub>2</sub> footprint



Made of 75% recycled aluminium



Discover our tools in action through the Interactive Product videos on our website: [www.lopital.com](http://www.lopital.com)

Art. nr.	Article description	Inside length	Inside width
6720 0005	Mattress Basic 186 cm	170	70
6720 0015	Mattress Basic 206 cm	190	70
6720 0010	Mattress Deluxe 186 cm	170	70
6720 0020	Mattress Deluxe 206 cm	190	70
6720 0030	Mattress Deluxe 226 cm	170	70



Suitable for people with severely reduced to very severely reduced mobility

**Flame retardant (M2 standard).**

Meets M2 fire safety standards for added protection in professional care environments.

**Soft Touch mattress.**


Extra comfortable, soft to the touch and equipped with a non-slip coating for added safety.

**High scratch resistance**

The durable surface resists wear and tear, maintaining a clean and professional appearance over time.

**Antibacterial.**

The mattress is antibacterial due to the addition of silver ions.

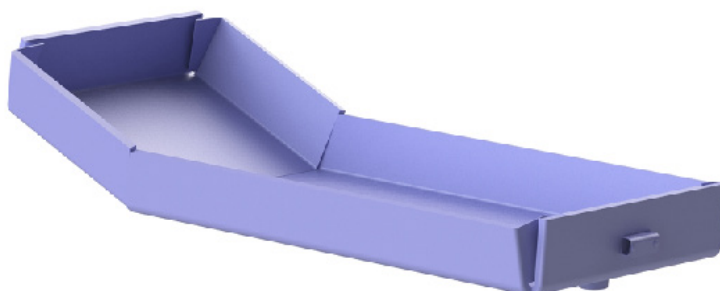**REACH-compliant & cadmium-free.**

Produced in accordance with REACH standards and free from cadmium for safer, environmentally responsible use.

**Safe & skin-friendly mattress.**

Free from harmful substances like phthalates – ideal for sensitive skin and safe daily use.

**Durable multi-layer mattress.**

Composed of four layers of PVC for enhanced comfort, long-lasting support and excellent resistance to water and wear.

