

Elexo XXL

Strength in care, comfort in every move

5100 5300

Elexo XXL shower toilet chair

The Elexo XXL is specifically engineered to meet the demands of clients with higher weight and size requirements, offering a stable, robust solution for hygienic care. With its wide, ergonomically shaped seat and extended, fully foldable armrests, it supports comfort and safe transfers. The swing-out, removable footrests and integrated calf support allow caregivers to position clients precisely, while minimizing awkward handling and physical strain. Its generous dimensions and thoughtful design make it a dependable choice in bariatric care settings.

Height adjustability and tilt capability elevate the functional flexibility of the Elexo XXL. The electric height range from low to high allows caregivers to maintain optimal ergonomic postures during use. At its highest position, the chair achieves a 20° tilt for improved seating balance and client support, facilitating both care and comfort. The dual-motor drive ensures smooth, stable adjustments, while integrated safety features like emergency lowering and sleep mode enhance reliability in daily operation.

Ergonomic Use. Equipped with smart handle and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.

Reliable & Long-Lasting. Built for daily care and backed by a 60-month warranty, offering a dependable performance.

Antibacterial & Corrosion-Resistant. The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

 Discover our tools in action through the Interactive Product videos on our website: www.lopital.com



Technical specifications

Max. load	220 kg / 485 lbs
Product weight	71 kg / 156.5 lbs
Seat width	70 cm / 27.6 inch
Seat depth	54 cm / 21.3 inch
Seat height (adjustable)	50 to 105 cm / 19.7 to 41.3 inch



Suitable for people with
very slightly to very
severely reduced mobility


Ergonomic handle
for smooth and light
operation

Guide rails for
bedpan or toilet
bucket

Strong and stable
frame for maximum
safety and support

Foldable and extended
armrests

Ergonomic seating comfort

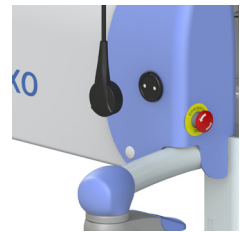

Charging point with
magnetic coupling,
no separate battery
required.



Independently foldable
armrests.

Bumper to prevent
damage to the frame



Optional

6960 8026
Comfort leg support
set.



6840 5500
Adjustable head/
neck support.



Options

- Product no. 6840 5500 Head/neck support
- Product no. 6950 1104 Bedpan with lid plastic
- Product no. 6975 8044 Seatbelt 2-point
- Product no. 6975 8045 Seatbelt 4-point
- Product no. 6960 8026 Comfort leg support set

Technical specifications

Total length base frame	119 cm / 46.9 inch
Total width base frame	89 cm / 35 inch
Turning radius	140 cm / 55.1 inch
Tilt setting	0 to 20 degrees
Electrical power	24V
Battery capacity	7Ah
Waterproofness	IPX4



For more information also visit our website: www.lopital.com
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk
Tel +31 (0)13 52.39.300, E-mail info@lopital.nl, www.lopital.nl



ISO CERTIFIED
9001 - 13485 - 14001

MADE IN
THE NETHERLANDS



Scan for
information

All of the aids produced by Lopital comply with the applicable standards for Class 1 medical products (MDR (EU) 2017/745).

Copyright © 2025 Lopital. All rights reserved.