

# Elexo XXL

## Strength in care, comfort in every move

5100 5300

### Elexo XXL shower toilet chair

The Elexo XXL is specifically engineered to meet the demands of clients with higher weight and size requirements, offering a stable, robust solution for hygienic care. With its wide, ergonomically shaped seat and extended, fully foldable armrests, it supports comfort and safe transfers. The swing-out, removable footrests and integrated calf support allow caregivers to position clients precisely, while minimizing awkward handling and physical strain. Its generous dimensions and thoughtful design make it a dependable choice in bariatric care settings.

Height adjustability and tilt capability elevate the functional flexibility of the Elexo XXL. The electric height range from low to high allows caregivers to maintain optimal ergonomic postures during use. At its highest position, the chair achieves a 20° tilt for improved seating balance and client support, facilitating both care and comfort. The dual-motor drive ensures smooth, stable adjustments, while integrated safety features like emergency lowering and sleep mode enhance reliability in daily operation.

**Ergonomic Use.** Equipped with smart handle and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.

**Reliable & Long-Lasting.** Built for daily care and backed by a 60-month warranty, offering a dependable performance.

**Antibacterial & Corrosion-Resistant.** The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

 Discover our tools in action through the Interactive Product videos on our website: [www.lopital.com](http://www.lopital.com)



#### Technical specifications

|                          |                                  |
|--------------------------|----------------------------------|
| Max. load                | 220 kg / 485 lbs                 |
| Product weight           | 71 kg / 156.5 lbs                |
| Seat width               | 70 cm / 27.6 inch                |
| Seat depth               | 54 cm / 21.3 inch                |
| Seat height (adjustable) | 50 to 105 cm / 19.7 to 41.3 inch |



70550007  
Hand control  
with battery  
indication



Suitable for people with  
very slightly to very  
severely reduced mobility


Ergonomic handle  
for smooth and light  
operation

Guide rails for  
bedpan or toilet  
bucket

Strong and stable  
frame for maximum  
safety and support

Foldable and extended  
armrests

Ergonomic seating comfort

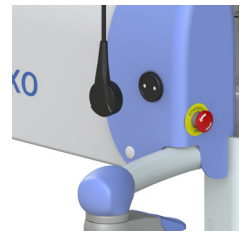

Charging point with  
magnetic coupling,  
no separate battery  
required.



Independently foldable  
armrests.

Bumper to prevent  
damage to the frame



## Optional

6960 8026  
Comfort leg support  
set.



6840 5500  
Adjustable head/  
neck support.



### Options

- 6840 5500 Head/neck support
- 6950 1104 Bedpan with lid plastic
- 6960 8026 Comfort leg support set

### Technical specifications

|                         |                    |
|-------------------------|--------------------|
| Total length base frame | 119 cm / 46.9 inch |
| Total width base frame  | 89 cm / 35 inch    |
| Turning radius          | 140 cm / 55.1 inch |
| Tilt setting            | 0 to 20 degrees    |
| Electrical power        | 24V                |
| Battery capacity        | 7Ah                |
| Waterproofness          | IPX4               |

For more information also visit our website: [www.lopital.com](http://www.lopital.com)  
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk  
Tel +31 (0)13 52.39.300, E-mail [info@lopital.nl](mailto:info@lopital.nl), [www.lopital.nl](http://www.lopital.nl)



ISO CERTIFIED  
9001 - 13485 - 14001

MADE IN  
THE NETHERLANDS



Scan for  
information

All of the aids produced by Lopital comply with the applicable standards for Class 1 medical products (MDR (EU) 2017/745).

Copyright © 2025 Lopital. All rights reserved.