

# Lotus

## Empowering safe transfers with comfort and control

5900 8200 / 8250 / 8270

### Lotus active standup aid

The Lotus is an active standup aid designed for individuals who can still stand independently and possess sufficient core stability but may have difficulty walking or are unable to walk. Caregivers can safely and efficiently move the client after they have seated themselves on the Lotus.

The Lotus is intended for short transfers, such as moving from the bed to the toilet. Because the client assumes an active seated position, less effort is required to return to a standing position. Optimal comfort is provided for both the client and the caregiver.

The Lotus features swivelling seat components equipped with snap lock closure, comfortable lower leg support, and a leg spread function. The leg spread function allows the Lotus to be brought close to the client, enabling them to actively take a seat. The Lotus encourages the client to actively participate in transfers, which helps train muscles and promotes independence.

**Ergonomic Use.** Equipped with smart handles and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.

**Reliable & Long-Lasting.** Built for daily care and backed by a 60-month warranty, offering a dependable performance.

**Antibacterial & Corrosion-Resistant.** The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

 Discover our tools in action through the Interactive Product videos on our website: [www.lopital.com](http://www.lopital.com)



Suitable for people with very slightly to moderately reduced mobility

#### Technical specifications

|                |                                  |
|----------------|----------------------------------|
| Max. load      | 200 - 227 kg / 440.9 - 500.4 lbs |
| Product weight | 32 kg / 70.5 lbs                 |
| Warranty       | 5 years                          |

5900 8200  
Lotus active standup aid

Ergonomic seat  
comfort

Lower leg support

Footplate

Pedal

Handle

Push handle


Swivelling seat  
components equipped  
with snap lock closure

Mechanical base  
spreading

Wheel cover

Swivel wheel Ø100

| Art.no.   | Art. description                    | Material | Lifting<br>capacity | Length x<br>width x<br>height                 | Width legs   | Mechanical base<br>spreading   | Weight          | Client height                |
|-----------|-------------------------------------|----------|---------------------|---|--------------|--------------------------------|-----------------|------------------------------|
| 5900 8200 | Lotus active standup aid            | steel    | 200 kg<br>440 lbs   | 96 cm x 66 cm x 110 cm<br>37.8" x 26" x 43.3" | 66 cm<br>26" | 49 cm - 86 cm<br>19.3" - 33.9" | 32 kg<br>70 lbs | 170 cm - 195 cm<br>67" - 77" |
| 5900 8250 | Lotus Compact active<br>standup aid | steel    | 200 kg<br>440 lbs   | 96 cm x 66 cm x 101 cm<br>37.8" x 26" x 39.8" | 66 cm<br>26" | 49 cm - 86 cm<br>19.3" - 33.9" | 32 kg<br>70 lbs | 135 cm - 170 cm<br>53" - 67" |
| 5900 8270 | Lotus XL active standup aid         | steel    | 227 kg<br>500 lbs   | 96 cm x 71 cm x 110 cm<br>37.8" x 28" x 43.3" | 66 cm<br>26" | 49 cm - 86 cm<br>19.3" - 33.9" | 32 kg<br>70 lbs | 170 cm - 195 cm<br>67" - 77" |

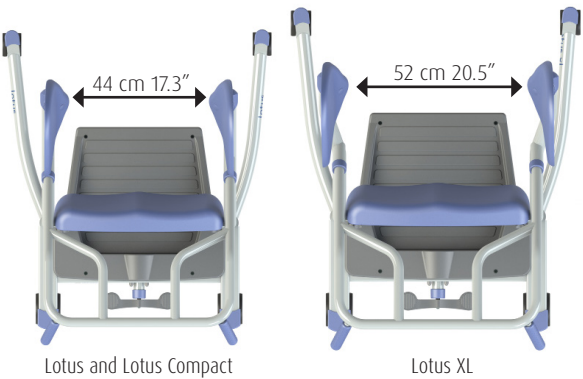

Lotus and Lotus Compact

Lotus XL

LFT8200-2025-01 EN



For more information also visit our website: [www.lopital.com](http://www.lopital.com)  
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk  
Tel +31 (0)13 52.39.300, E-mail [info@lopital.nl](mailto:info@lopital.nl), [www.lopital.nl](http://www.lopital.nl)

  
**Green Deal**

 **ISO CERTIFIED**  
9001 - 13485 - 14001

 **MADE IN  
THE NETHERLANDS**



  
Scan for  
information

All of the aids produced by Lopital comply with the applicable standards for Class 1 medical products (MDR (EU) 2017/745).

Copyright ©2025 Lopital. All rights reserved.