

Lotus

Empowering safe transfers with comfort and control

5900 8200 / 8250 / 8270

Lotus active standup aid

The Lotus is an active standup aid designed for individuals who can still stand independently and possess sufficient core stability but may have difficulty walking or are unable to walk. Caregivers can safely and efficiently move the client after they have seated themselves on the Lotus.

The Lotus is intended for short transfers, such as moving from the bed to the toilet. Because the client assumes an active seated position, less effort is required to return to a standing position. Optimal comfort is provided for both the client and the caregiver.

The Lotus features swivelling seat components equipped with snap lock closure, comfortable lower leg support, and a leg spread function. The leg spread function allows the Lotus to be brought close to the client, enabling them to actively take a seat. The Lotus encourages the client to actively participate in transfers, which helps train muscles and promotes independence.

Ergonomic Use. Equipped with smart handles and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.

Reliable & Long-Lasting. Built for daily care and backed by a 60-month warranty, offering a dependable performance.

Antibacterial & Corrosion-Resistant. The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

 Discover our tools in action through the Interactive Product videos on our website: www.lopital.com



Technical specifications

Max. load	200 - 227 kg / 440.9 - 500.4 lbs
Product weight	32 - 36 kg / 70.5 - 79.4 lbs
Warranty	5 years



Suitable for people with very slightly to moderately reduced mobility

5900 8200

Lotus active standup aid

Ergonomic seat comfort

Handle

Push handle

Lower leg support


Footplate

Mechanical base spreading

Wheel cover


Pedal

Swivel wheel Ø100

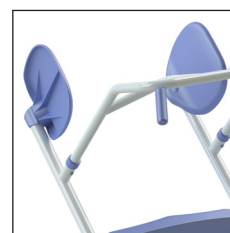

Swivelling seat components equipped with snap lock closure

Art.no.	Art. description	Material	Lifting capacity	Length x width x height	Width legs	Mechanical base spreading	Weight	Client height
5900 8200	Lotus active standup aid	steel	200 kg 440 lbs	96 cm x 66 cm x 110 cm 37.8" x 26" x 43.3"	66 cm 26"	49 cm – 86 cm 19.3" – 33.9"	32 kg 70 lbs	170 cm – 195 cm 67" – 77"
5900 8250	Lotus Compact active standup aid	steel	200 kg 440 lbs	96 cm x 66 cm x 101 cm 37.8" x 26" x 39.8"	66 cm 26"	49 cm – 86 cm 19.3" – 33.9"	32 kg 70 lbs	135 cm – 170 cm 53" – 67"
5900 8270	Lotus XL active standup aid	steel	227 kg 500 lbs	102 cm x 74 cm x 110 cm 40.2" x 29.1" x 43.3"	74 cm 29.1"	57 cm – 94 cm 22.4" – 37"	36 kg 79.4 lbs	170 cm – 195 cm 67" – 77"

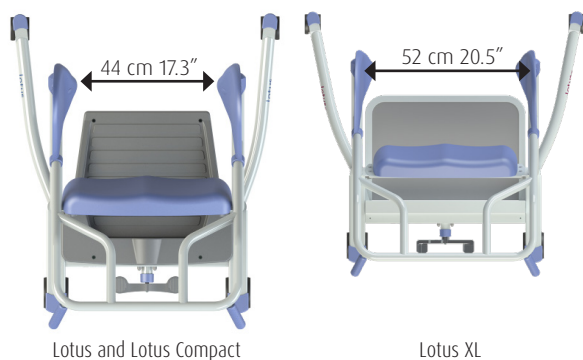

Lotus and Lotus Compact

Lotus XL

LFT8200-2025-01 EN



For more information also visit our website: www.lopital.com
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk
Tel +31 (0)13 52.39.300, E-mail info@lopital.nl, www.lopital.nl



ISO CERTIFIED
9001 - 13485 - 14001

MADE IN
THE NETHERLANDS



All of the aids produced by Lopital comply with the applicable standards for Class 1 medical products (MDR (EU) 2017/745).

Copyright ©2026 Lopital. All rights reserved.