

Reflex XXL

Empowering care, without limits

5100 5900

Reflex XXL shower toilet chair

The Reflex XXL extends the proven Reflex design philosophy to meet higher-capacity needs — now enhanced with adjustable ergonomics specifically developed for bariatric care. The height-adjustable backrest and width-adjustable armrests provide personalised support for a wide range of body shapes. The armrests lock securely and can be released effortlessly with a push button, enabling smooth and controlled transfers. With its robust stainless-steel frame, reinforced stability and powerful electric operation, the Reflex XXL delivers uncompromised comfort, safety and confidence for larger users.

Purpose-built for demanding care settings, the Reflex XXL combines reliability, comfort and advanced functionality in one durable solution. Its smooth, synchronised height and tilt adjustments ensure precise positioning, including a forward negative-tilt function that creates an intuitive stand-up assist. Together with its spacious design and high load capacity, the Reflex XXL supports safe, dignified and efficient care — for both caregiver and client.

Adaptive Support. The height-adjustable backrest and width-adjustable armrests provide optimal support for a wide range of body profiles. The armrests lock securely into place and can be released with a single push button, ensuring smooth and controlled transfers.

Reliable & Long-Lasting. Built for daily care and backed by a 60-month warranty, offering a dependable performance.

Antibacterial & Corrosion-resistant. The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

Assistive Tilt Function. Thanks to the negative seat tilt, the Reflex XXL offers an effective stand-up support function. This enhances user independence while reducing physical strain for the caregiver.




 Discover our tools in action through the Interactive Product videos on our website: www.lopital.com

Technical specifications

Max. load	275 kg / 606.3 lbs
Product weight	70 kg / 154.3 lbs
Seat width (adjustable)	70 – 80 cm / 27.6 – 31.5 inch
Seat depth	53 cm / 20.9 inch
Seat height (adjustable)	45 – 76 cm / 19.7 – 29.9 inch



Head/neck support.

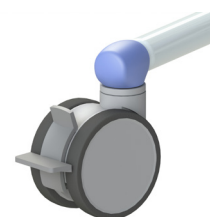

3 Swivel castors and 1 directional castor with individual brakes.



Suitable for people with moderately to severely reduced mobility

5100 5900

Reflex XXL shower toilet chair

Optional

6950 1104 Bedpan and lid plastic.



Technical specifications

Total length base frame	121 cm / 47.6 inch
Total width base frame	84 cm / 33.1 inch
Turning radius	140 cm / 55.1 inch
Tilt setting	0 to 30 degrees backward, 10 degrees forward
Electrical power	24V
Battery capacity	7Ah
Charger (Mascot 3743)	110-240V AC, 50-60Hz, max. 0.35A, cable length 1.5 m
Head/neck support depth/height	depth 7 cm / 2.8 inch; height 12 cm / 4.7 inch
Footrest height adjustment	8 cm / 3.2 inch

LFT5900-2026-01 EN



For more information also visit our website: www.lopital.com
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk
Tel +31 (0)13 52.39.300, E-mail info@lopital.nl, www.lopital.nl



ISO CERTIFIED
9001 - 13485 - 14001

MADE IN
THE NETHERLANDS



All of the aids produced by Lopital comply with the applicable standards for Class I medical products (MDR (EU) 2017/745).

Copyright ©2026 Lopital. All rights reserved.