

Tango

Comfort and flexibility in daily care

5100 5700

Tango shower toilet chair

The Tango is designed for everyday care moments where comfort, safety and efficiency are key. Its stable frame combined with a compact, manoeuvrable design offers practical confidence during showers, toilet visits and transfers – with a focus on both safety and comfort.

The ergonomic design with smart handles makes the chair light and manoeuvrable to use, significantly reducing physical strain for the caregiver. For the client, the folding and swing-away footrests ensure comfortable and safe entry and exit.

The four braked castors ensure easy and stable movement on various floor types. Thanks to the use of high-quality, easy-care materials, the chair is also easy to clean - hygiene and durability guaranteed.

Lopital Tango – moving with care, wherever care moves.

Compact & Manoeuvrable Design. With its lightweight and compact frame, Tango is easy to handle in smaller spaces - ideal for home care and tighter bathrooms.

Ergonomic Use. Equipped with smart handles and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.

Reliable & Long-Lasting. Built for daily care and backed by a 60-month warranty, offering a dependable performance.

Antibacterial & Corrosion-Resistant. The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

 Discover our tools in action through the Interactive Product videos on our website: www.lopital.com



Technical specifications

Max. load	180 kg / 400 lbs
Product weight	19 kg / 41.9 lbs
Seat width	47 cm / 18.5 inch
Seat depth	44 cm / 17.3 inch



Suitable for people with very slightly to severely reduced mobility




Optional

6950 1104
Plastic bedpan
and lid.



6975 8044
Safety belt 2-point.

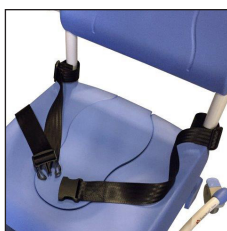

6840 5660
Front locking
armrests.



Options

- Product no. 6840 5660 Armrests with front locking set
- Product no. 6840 5661 Armrests with front dual locking set
- Product no. 6950 1104 Bedpan with lid plastic
- Product no. 6975 8044 Seatbelt 2-point
- Product no. 6975 8045 Seatbelt 4-point
- Product no. 6950 5624 Platilon backrest
- Product no. 6950 5625 Platilon seat U-model
- Product no. 6960 8026 Comfort leg support set
- Product no. 7500 5693 Raising set 10 up to 30 mm
- Product no. 7057 0508 Castor 125 mm Directional

Technical specifications

Seat height	55 cm / 21.7 inch
Total length base frame	84 cm / 33.1 inch
Total width base frame	66 cm / 26.0 inch
Turning radius	98 cm / 38.6 inch

LFT5700-2025-01 EN



For more information also visit our website: www.lopital.com
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk
Tel +31 (0)13 52.39.300, E-mail info@lopital.nl, www.lopital.nl

