

# Elexo

Comfort, control and dignity – in every care movement

5100 5500

## Elexo shower toilet chair

The Elexo is a fully electric height-adjustable shower toilet chair designed to support clients who require comfort, safety, and dignity during personal care routines. With its quiet motorised functions, the Elexo allows for seamless positioning transitions – enabling caregivers to maintain safe and efficient workflows without disrupting the client's sense of stability.

The chair is suited for use in both shower and toilet settings and can be easily positioned over most toilet models. Its compact design allows for smooth navigation in tight spaces, making it ideal for modern care environments. Integrated user-friendly controls and thoughtful design elements ensure that both client and caregiver benefit from a smooth, secure care experience.

The Elexo is a trusted solution for professionals seeking to elevate care standards with precision, ease, and adaptability.

**Ergonomic Use.** Equipped with smart handles and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.


**Reliable & Long-Lasting.** Built for daily care and backed by a 60-month warranty, offering a dependable performance.

**Antibacterial & Corrosion-Resistant.** The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

 Discover our tools in action through the Interactive Product videos on our website: [www.lopital.com](http://www.lopital.com)

### Technical specifications

Max. load	180 kg / 400 lbs
Product weight	47 kg / 103.6 lbs
Seat width	47 cm / 18.5 inch
Seat depth	44 cm / 17.3 inch
Seat height (adjustable)	50 to 105 cm / 19.7 to 41.3 inch

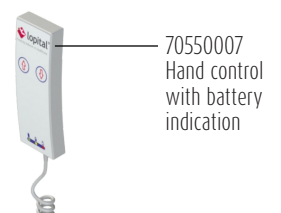

Suitable for people with  
very slightly to very  
severely reduced mobility

Ergonomic handle for smooth and light operation

Armrests with front closure for safety and security



Charging point with magnetic coupling, no separate battery required.

Guide rails for bedpan or toilet bucket

Ergonomic seating comfort

Strong and stable frame for maximum safety and support

Bumper to prevent damage to the frame



## Optional

6960 8026  
Comfort leg support set.



6840 5661  
Dual locking armrests.



6840 5500  
Adjustable head/neck support.



### Options

- Product no. 6830 5650 Straight armrests
- Product no. 6840 5661 Armrests with front dual locking set
- Product no. 6950 1104 Bedpan with lid plastic
- Product no. 6975 8044 Seatbelt 2-point
- Product no. 6975 8045 Seatbelt 4-point
- Product no. 6950 5624 Platilon backrest
- Product no. 6950 5625 Platilon seat U-model
- Product no. 6960 8026 Comfort leg support set
- Product no. 7500 5693 Raising set 10 up to 30 mm
- Product no. 7057 0508 Castor 125 mm Directional

### Technical specifications

Total length base frame	104 cm / 40.9 inch
Total width base frame	70 cm / 27.6 inch
Turning radius	118 cm / 46.5 inch
Tilt setting	0 to 20 degrees
Electrical power	24V
Battery capacity	7Ah
Waterproofness	IPX4

For more information also visit our website: [www.lopital.com](http://www.lopital.com)  
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk  
Tel +31 (0)13 52.39.300, E-mail [info@lopital.nl](mailto:info@lopital.nl), [www.lopital.nl](http://www.lopital.nl)



ISO CERTIFIED  
9001 - 13485 - 14001

MADE IN  
THE NETHERLANDS



Scan for  
information