

# Nemo

Safe and secure tilting comfort for any care situation

5100 5100

## Nemo shower toilet chair

The Nemo is a castored and tiltable shower and toilet chair specially designed for clients who need extra support during daily care. The manually tilting shower and toilet chair allows clients to be placed in a comfortable and stable position, while care professionals continue to work in an ergonomically sound position. This contributes to physical relief and increases safety for both parties.

The chair is suitable for use both in the shower and as a toilet seat and can be easily placed over the toilet. With its compact size and central braking system, the Nemo is highly manoeuvrable and safe to use, even in smaller spaces. The sturdy, low-maintenance construction is designed for daily intensive use.

The Nemo is a reliable aid that contributes to high-quality daily care with attention and ease of use.

**Manually tiltable** up to 30 degrees by gas spring for optimal positioning and ergonomic care.

**Ergonomic Use.** Equipped with smart handles and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.

**Reliable & Long-Lasting.** Built for daily care and backed by a 60-month warranty, offering a dependable performance.

**Antibacterial & Corrosion-Resistant.** The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

### Options

- Product no. 6840 5660 Armrests with front locking set
- Product no. 6840 5661 Armrests with front dual locking set
- Product no. 6950 1104 Bedpan with lid plastic
- Product no. 6975 8044 Seatbelt 2-point
- Product no. 6975 8045 Seatbelt 4-point
- Product no. 6950 5624 Platilon backrest
- Product no. 6950 5625 Platilon seat U-model
- Product no. 6960 8026 Comfort leg support set
- Product no. 7500 5693 Raising set 10 up to 30 mm
- Product no. 7057 0508 Castor 125 mm Directional



### Technical specifications

Max. load	150 kg / 330 lbs
Product weight	24 kg / 52.9 lbs
Seat width	47 cm / 18.5 inch
Seat depth	44 cm / 17.3 inch
Seat height	54 cm / 21.3 inch
Total length base frame	73 cm / 28.7 inch



Suitable for people with very slightly to severely reduced mobility





Discover our tools in action through the Interactive Product videos on our website: [www.lopital.com](http://www.lopital.com)

Adjustable head/  
neck support

Ergonomic  
seating  
comfort

Gas spring for  
smooth tilt  
adjustment

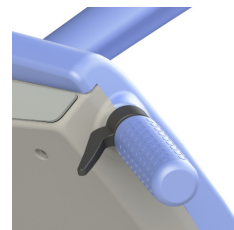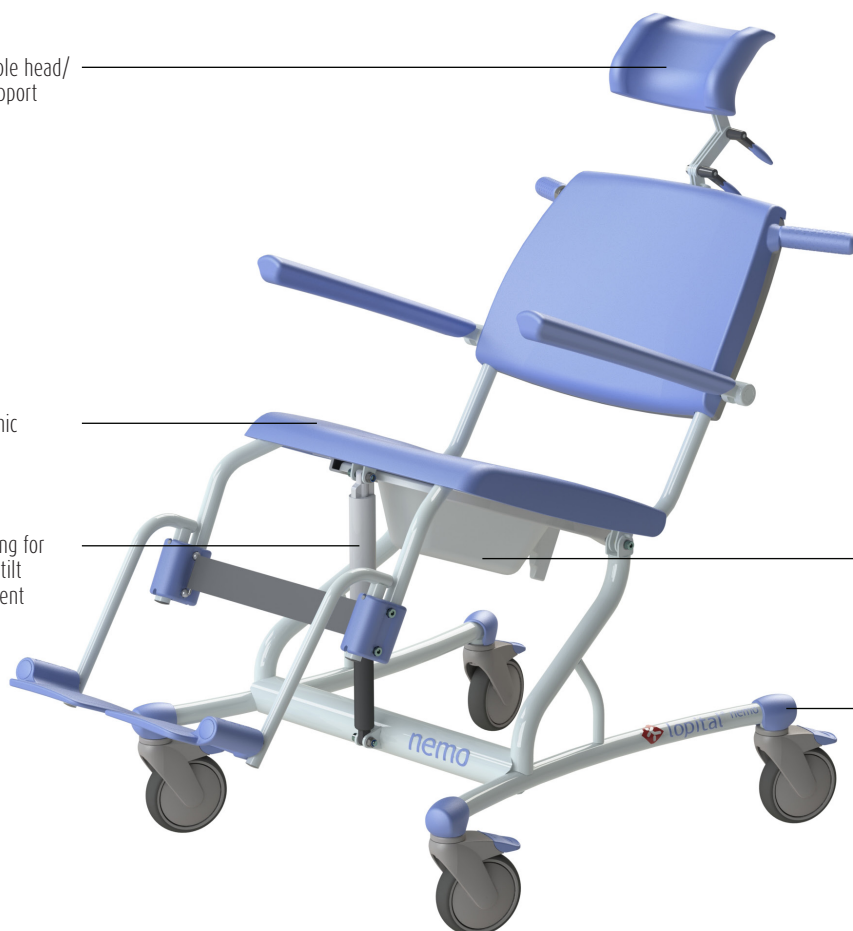

Gas spring handle for  
tilting position.

Guide rails for bedpan  
or toilet bucket

Bumper to prevent  
damage to the frame

## Optional


6960 8026  
Comfort leg support  
set.



6950 1104  
Bedpan with lid  
plastic.



6840 5500  
Adjustable head/neck  
support.

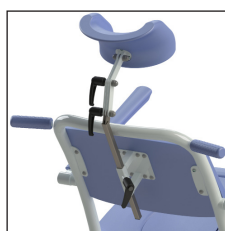

6840 5660  
Armrests with front  
locking set.



### Technical specifications

Total width basic frame	66 cm / 26 inch
Head / neck support depth / height	7 cm / 2.8 inch & 12 cm / 4.7 inch
Turning radius	112 cm / 44.1 inch
Tilt setting	0 to 30 degrees

LFT5100-2025-01 EN



For more information also visit our website: [www.lopital.com](http://www.lopital.com)  
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk  
Tel +31 (0)13 52.39.300, E-mail [info@lopital.nl](mailto:info@lopital.nl), [www.lopital.nl](http://www.lopital.nl)



ISO CERTIFIED  
9001 - 13485 - 14001

MADE IN  
THE NETHERLANDS



Scan for  
information

All of the aids produced by Lopital comply with the applicable standards for Class 1 medical products (MDR (EU) 2017/745).

Copyright © 2025 Lopital. All rights reserved.