

Nemo

Safe and secure tilting comfort for any care situation

5100 5100

Nemo shower toilet chair

The Nemo is a castored and tiltable shower and toilet chair specially designed for clients who need extra support during daily care. The manually tilting shower and toilet chair allows clients to be placed in a comfortable and stable position, while care professionals continue to work in an ergonomically sound position. This contributes to physical relief and increases safety for both parties.

The chair is suitable for use both in the shower and as a toilet seat and can be easily placed over the toilet. With its compact size and central braking system, the Nemo is highly manoeuvrable and safe to use, even in smaller spaces. The sturdy, low-maintenance construction is designed for daily intensive use.

The Nemo is a reliable aid that contributes to high-quality daily care with attention and ease of use.

Manually tiltable up to 30 degrees by gas spring for optimal positioning and ergonomic care.

Ergonomic Use. Equipped with smart handles and four braked castors to simplify manoeuvrability, ensuring high stability on various surfaces and minimizing the strain on the caregiver's joints.

Reliable & Long-Lasting. Built for daily care and backed by a 60-month warranty, offering a dependable performance.

Antibacterial & Corrosion-Resistant. The antibacterial coating with silver-ion supports hygiene and infection control, while the corrosion-resistant frame with smooth surfaces ensures durability and easy cleaning.

Options

- 6840 5660 Armrests with front locking set
- 6840 5661 Armrests with front dual locking set
- 6950 1104 Bedpan with lid plastic
- 6975 8044 Seatbelt 2-point
- 6975 8045 Seatbelt 4-point
- 6950 5624 Plaiton backrest
- 6950 5625 Plaiton seat U-model
- 6960 8026 Comfort leg support set
- 7500 5693 Raising set 10 up to 30 mm
- 7057 0508 Castor 125 mm Directional



Technical specifications

Max. load	150 kg / 330 lbs
Product weight	24 kg / 52.9 lbs
Seat width	47 cm / 18.5 inch
Seat depth	44 cm / 17.3 inch
Seat height	54 cm / 21.3 inch
Total length base frame	73 cm / 28.7 inch



Suitable for people with very slightly to severely reduced mobility




Discover our tools in action through the Interactive Product videos on our website: www.lopital.com

Adjustable head/
neck support

Ergonomic
seating
comfort

Gas spring for
smooth tilt
adjustment

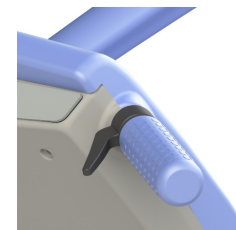

Gas spring handle for
tilting position.

Guide rails for bedpan
or toilet bucket


Bumper to prevent
damage to the frame

Optional


6960 8026
Comfort leg support
set.



6950 1104
Bedpan with lid
plastic.

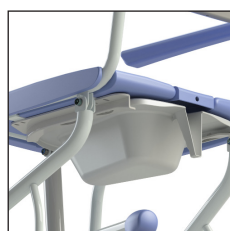

6840 5500
Adjustable head/neck
support.

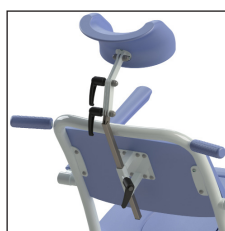

6840 5660
Armrests with front
locking set.



Technical specifications

Total width basic frame	66 cm / 26 inch
Head / neck support depth / height	7 cm / 2.8 inch & 12 cm / 4.7 inch
Turning radius	112 cm / 44.1 inch
Tilt setting	0 to 30 degrees

LFT5100-2025-01 EN



For more information also visit our website: www.lopital.com
Lopital Nederland B.V., Laarakkerweg 9, 5061 JR Oisterwijk
Tel +31 (0)13 52.39.300, E-mail info@lopital.nl, www.lopital.nl



ISO CERTIFIED
9001 - 13485 - 14001

MADE IN
THE NETHERLANDS



Scan for
information

All of the aids produced by Lopital comply with the applicable standards for Class 1 medical products (MDR (EU) 2017/745).

Copyright © 2025 Lopital. All rights reserved.